





ASAHI PCI Guide Wires

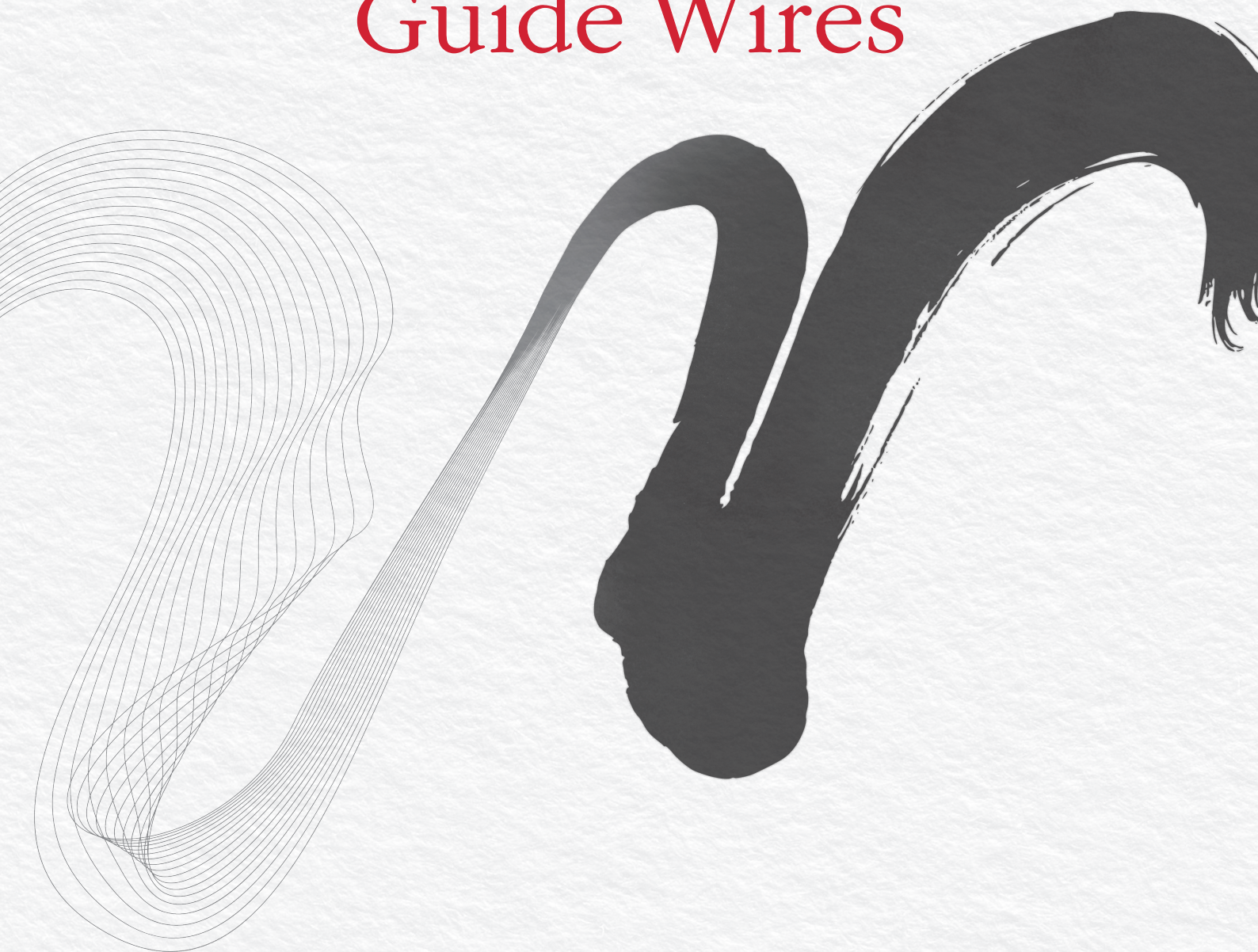

YOUR DREAMS. WOVEN TOGETHER.

Technology overcoming limitations


ASAHI INTECC's 4 core proprietary technologies incorporate multiple wire components to enhance guide wire durability and torque transmission.

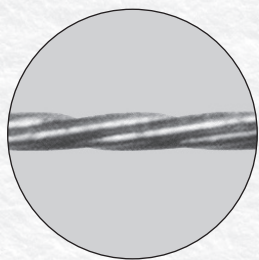

Twist Wire*


Provides tip durability and flexibility.

*Only for the following guide wires: SION blue, SION, SION black, SION blue ES, Gladius EX, Gladius MG, RG3, SUOH 03.

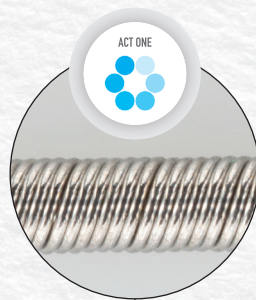



ACT ONE

ACT ONE is a multi-wire coil using our core technologies for wire drawing, wire forming, and torque. It optimizes the balance of torque force, torque response, durability, and rigidity in guidewires for high performance.

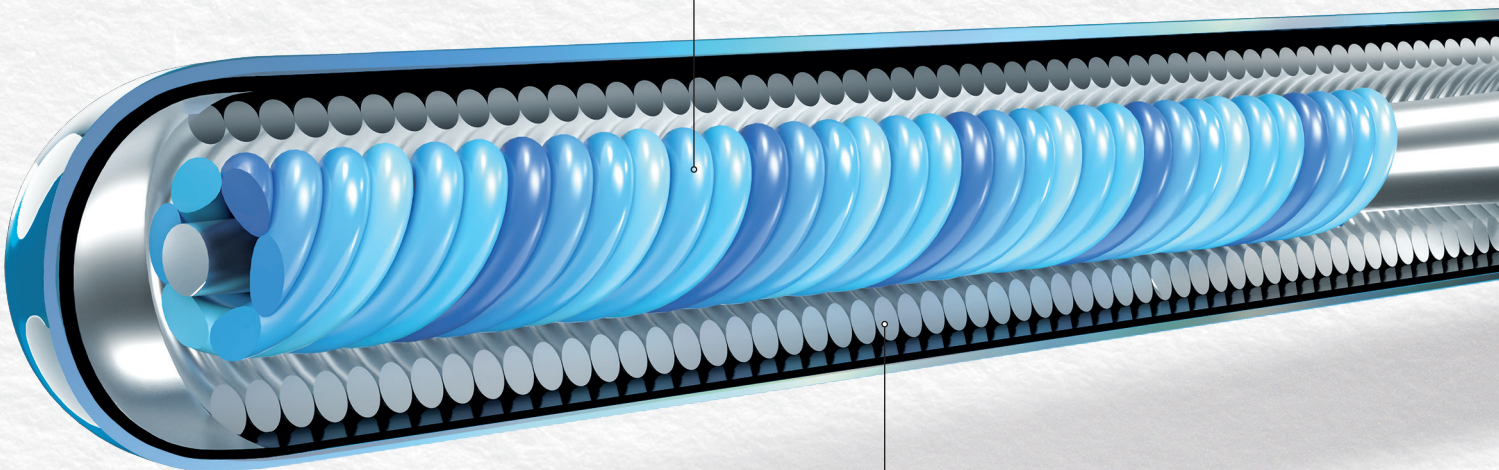

Tip Design

Variety of tip designs and profile to facilitate branch selection or to optimize penetration efficacy.

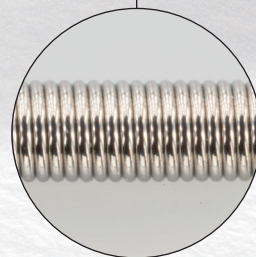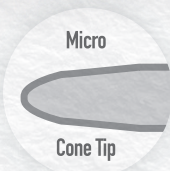


Joint-less spring coil*

Ensures a smooth transition from platinum to stainless steel spring coil and maintains the torqueability of the wire.



*Only for the following guide wires: SION blue, SION, SION black, SION blue ES, Gladius EX, Gladius MG, RG3, SUOH 03

ons

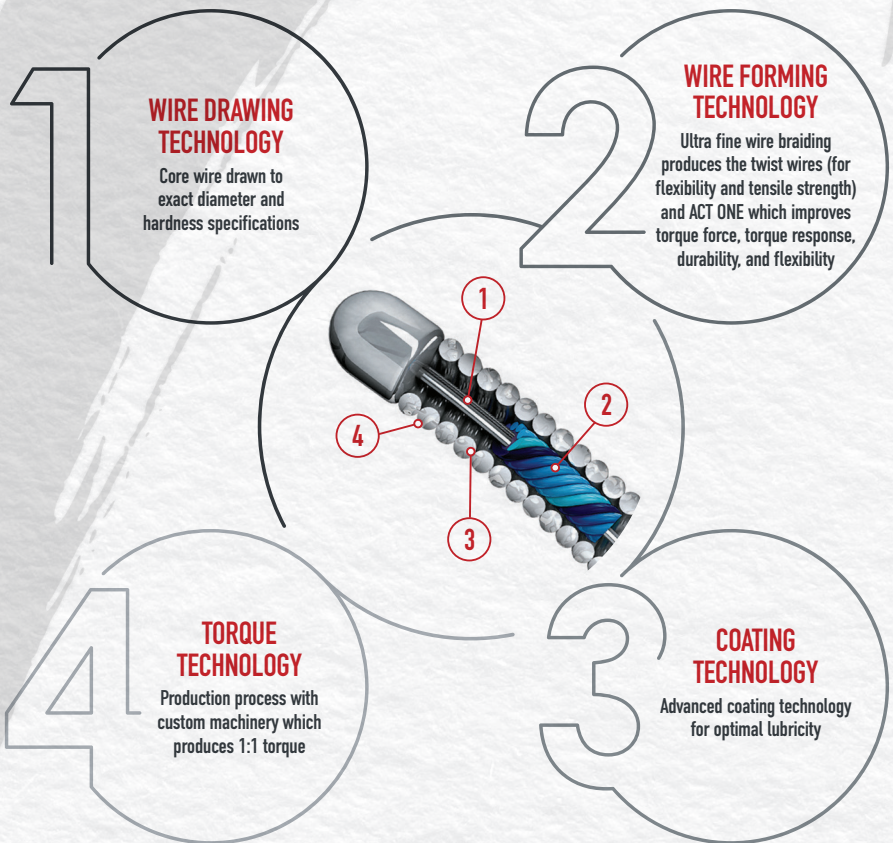


Stainless Steel Core

- Shaft support required for device delivery.
- Torque force required for operation.

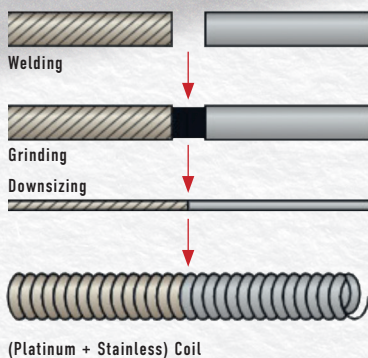
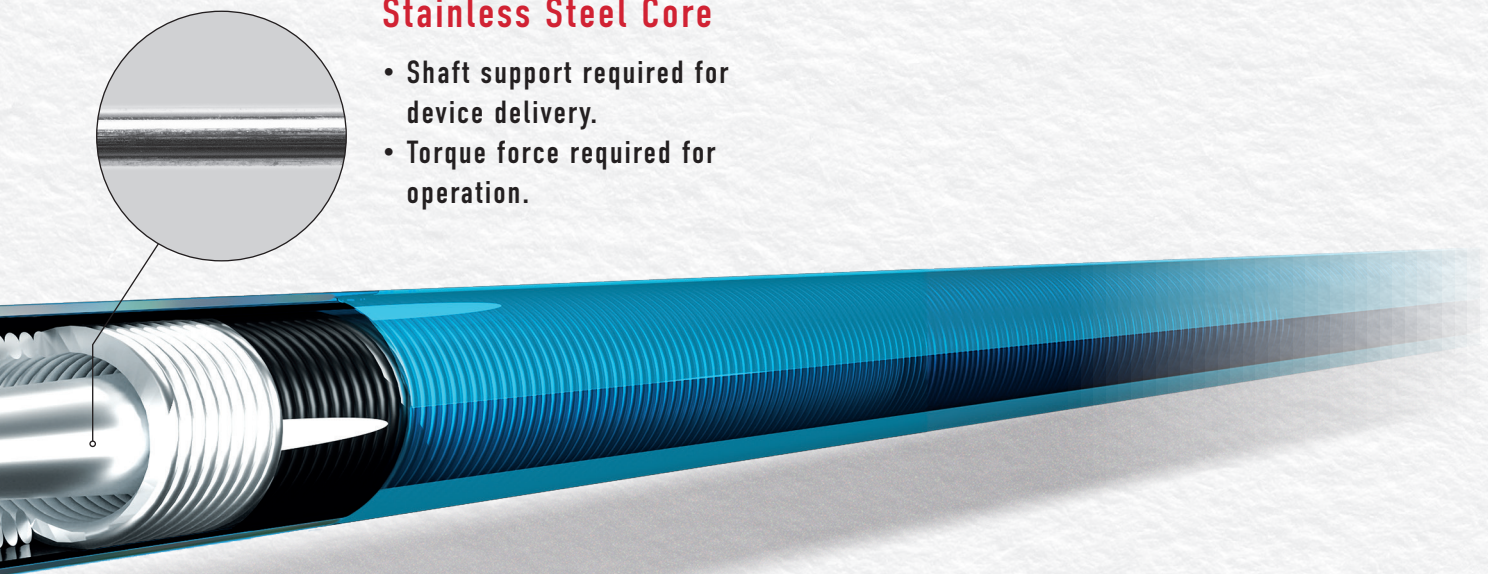


ASAHI Coating Technology

Variety of coating combinations to address any lesion type.


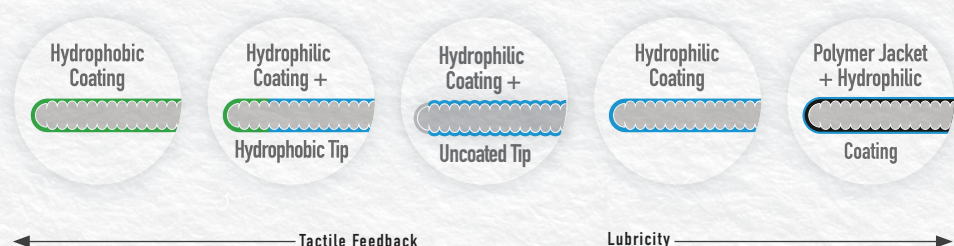


ASAHI SION blue

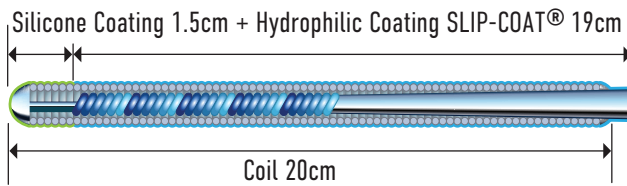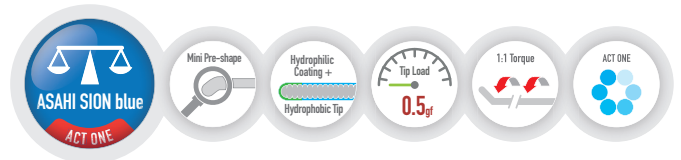



- First choice frontline guide wire.
- Balanced design which combines advanced flexibility and supportability.
- Silicone coating up to 15mm from the tip for safer procedures.

ASAHI SION

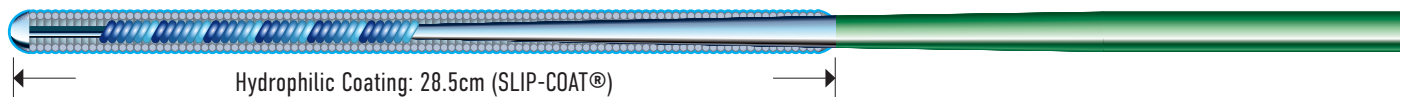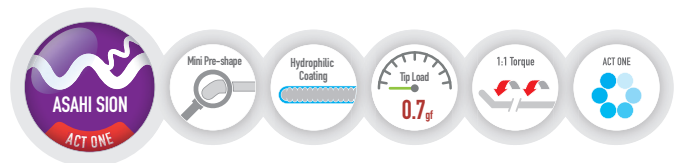



- Hydrophilic frontline guide wire.
- Flexible shaft and atraumatic tip for stress-free operation.

ASAHI SION black

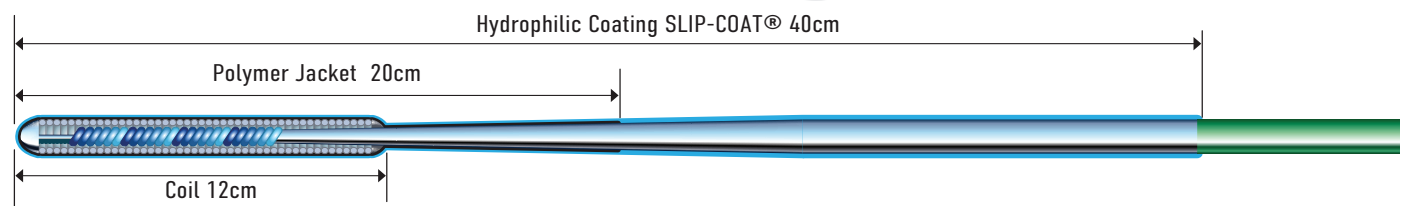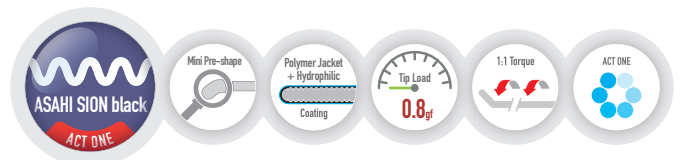


- Polymer jacket frontline guide wire.
- Designed to retain flexibility while crossing resistant stenosis and tortuous vessels.

ASAHI SION blue ES

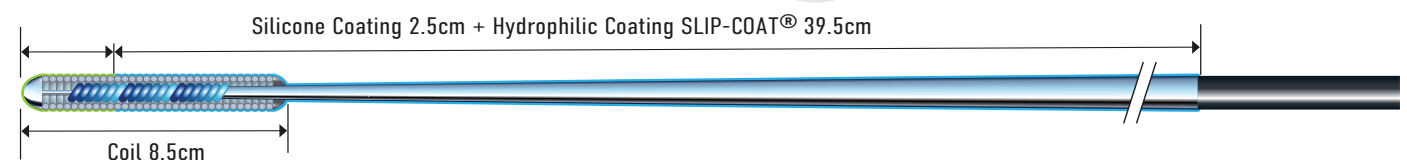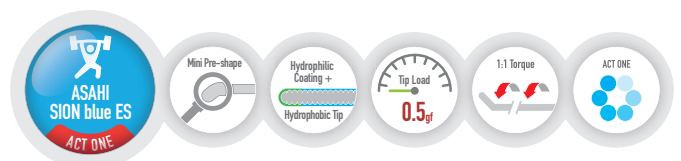

- Extra support frontline guide wire.
- Designed to stabilize the guide catheter for easier device delivery.
- Silicone coating up to 25mm from the tip for safer procedures.


Grand Slam

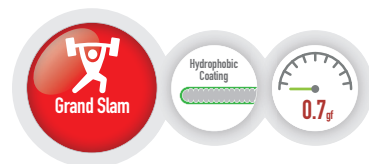

- Hydrophobic supportive guide wire.
- The stiff shaft provides extra support for device delivery in extreme cases.


Complex & Chronic Occlusion

Fielder XT-R

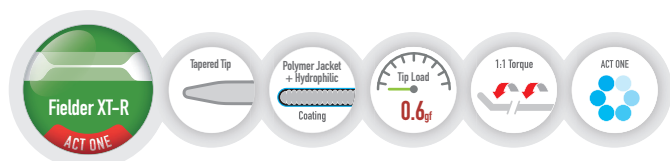



- Polymer jacket, tapered guide wire.
- First choice for CTO: recommended for loose tissue and narrow channel tracking with its low profile, flexible tip and lubricity.

Fielder XT-A

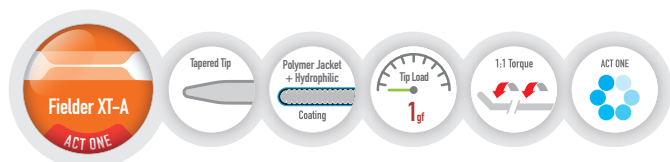



- Polymer jacket, tapered guide wire.
- First choice for CTO: low tip profile and lubricity for easier entry into the chronic occluded lesion.

ASAHI Miracle Neo 3

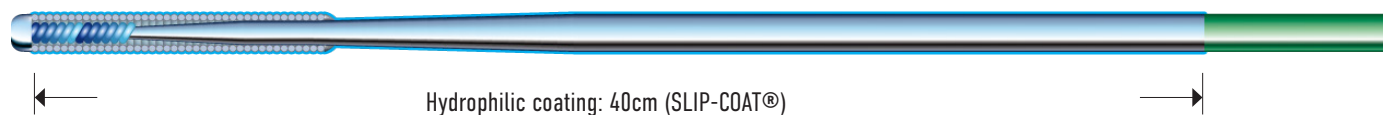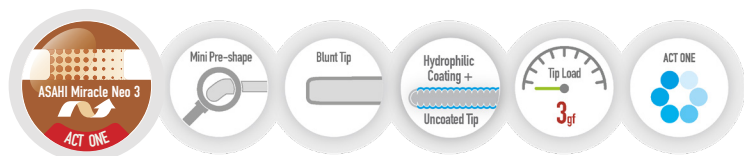



Complex & Chronic Occlusion

ASAHI Gladius MG

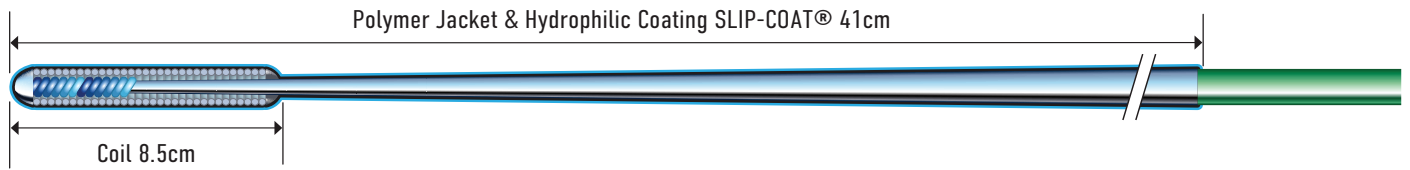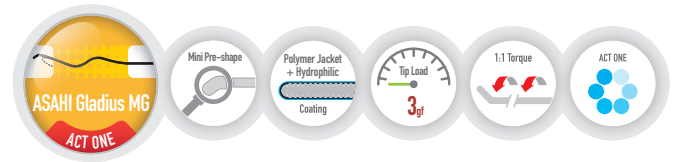


- Middleweight, non-tapered polymer jacket guide wire with modified distal core for narrow loop techniques and excellent lubricity for smooth navigation through long tortuous segments and ambiguous vessel course occlusions.

ASAHI Gladius EX

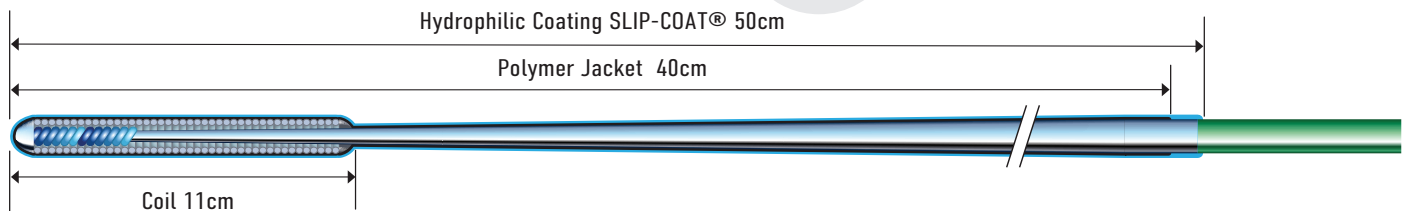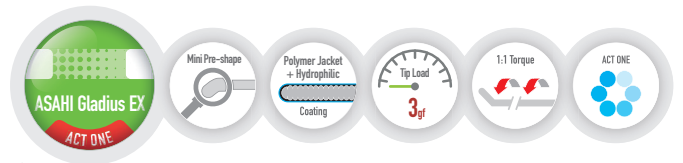

- Middleweight, non-tapered polymer jacket guide wire with excellent lubricity for smooth navigation through long tortuous segments and ambiguous vessel course occlusions.

ASAHI Gaia Next Family (Gaia Next 1, 2, 3)

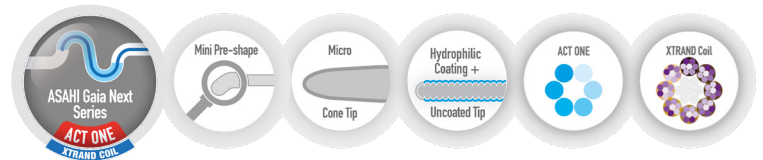



Tip Diameter:

ASAHI Gaia Next 1: 0.27mm (0.011 inch)

ASAHI Gaia Next 2: 0.30mm (0.012 inch)

ASAHI Gaia Next 3: 0.30mm (0.012 inch)

ASAHI Conquest Pro Series/ Confianza Pro Series

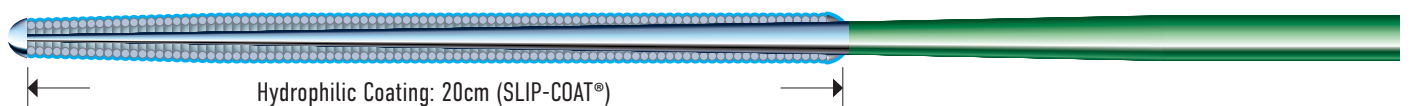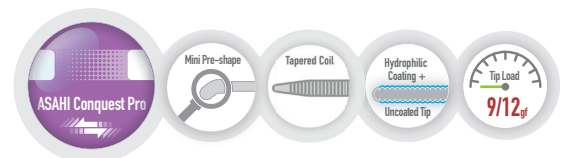

Tip Load:

ASAHI Conquest Pro / Confianza Pro: 9gf


ASAHI Conquest Pro 12 / Confianza Pro 12: 12gf

Tip Diameter:

ASAHI Conquest Pro / Confianza Pro: 0.23mm (0.009 inch)


ASAHI Conquest Pro 12 / Confianza Pro 12: 0.23mm (0.009 inch)

ASAHI Miracle 12/ MIRACLEbros 12

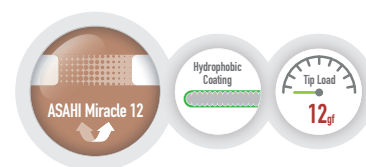

Externalization


ASAHI RG3

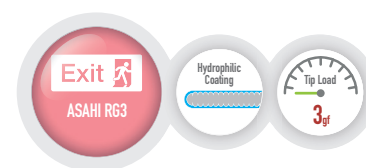

- Externalization guide wire with hydrophilic coating and 0.26mm (0.010 inch) shaft diameter which provides optimal flexibility, deliverability and stability for device advancement.



Retrograde

ASAHI SUOH 03

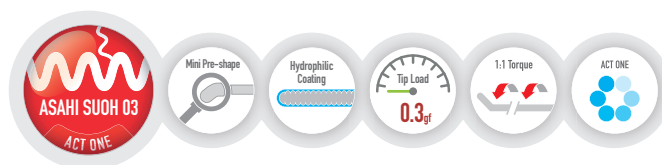


- Flexible tip provides superb trackability and safety while crossing extremely tortuous and angulated collaterals.

ASAHI Guide Wires Overview

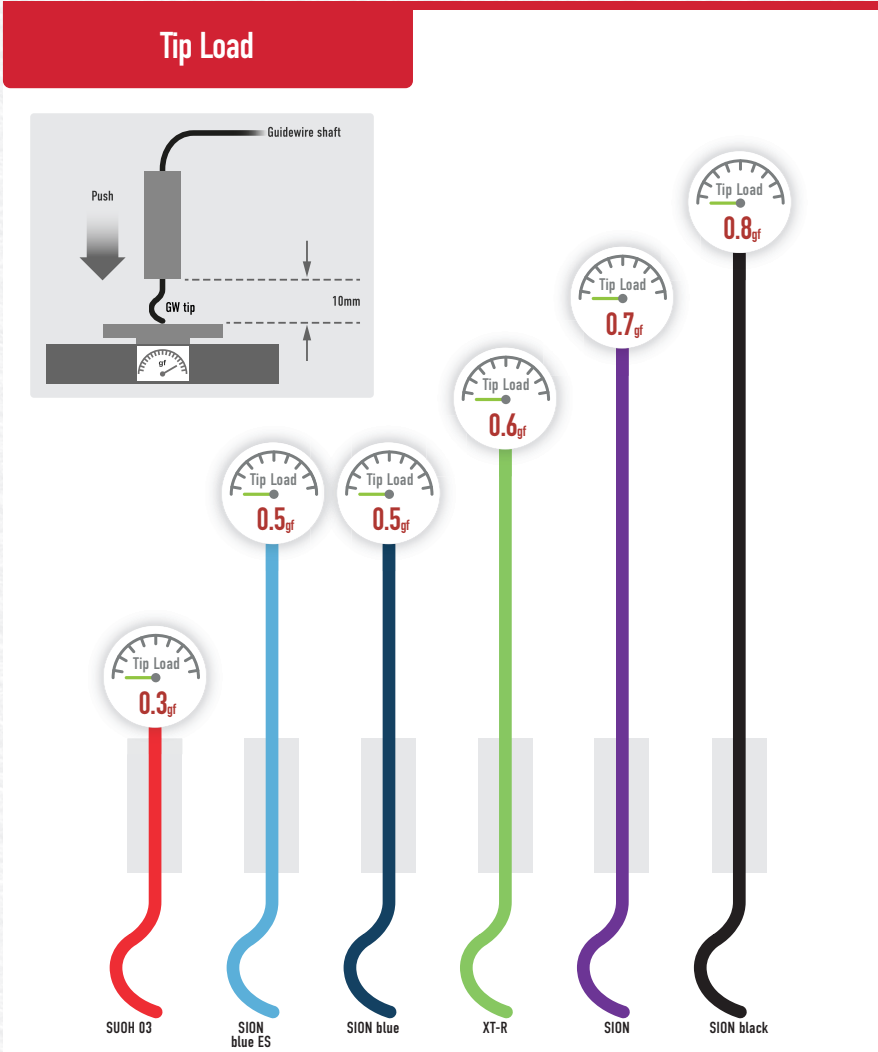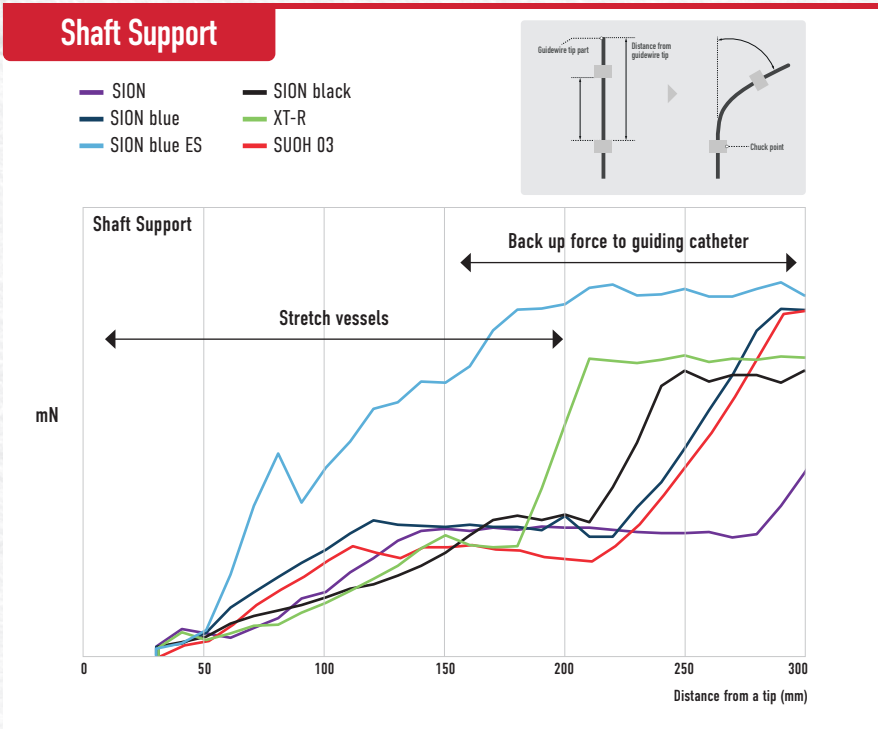


*All the data were obtained by company standardized test, which may differ from industry standardized tests.
 *All the data do not guarantee that all devices have exactly the same performance with the samples used for these tests.
 * The tip load values of our guide wires are the values assuming the force that the tip of the guide wire pushes the lesion on the basis of company standardized test. There are individual differences in values of actual devices.

Ordering Information & Lineup — Frontline Guide Wires

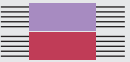Product	Catalog No.	Diameter (mm (inch))	Radiopaque Segment (cm)	Spring Coil Length (cm)	Length (cm)	Tip Shape	Clip Colour
ASAHI SION blue	AHW14R104S	0.36 (0.014)	3	20	190	Straight	
	AHW14R304S				300	Straight	
	AHW14R104J				190	J	
	AHW14R304J				300	J	
	AHW14R104P				190	Pre-shape	
	AHW14R304P				300	Pre-shape	
ASAHI SION	AHW14R101S	0.36 (0.014)	3	28	190	Straight	
	AHW14R301S				300	Straight	
	AHW14R101J				190	J	
	AHW14R301J				300	J	
	AHW14R101P				190	Pre-shape	
	AHW14R301P				300	Pre-shape	
ASAHI SION black	APW14R010S	0.36 (0.014)	3	12	190	Straight	
	APW14R310S				300	Straight	
	APW14R010J				190	J	
	APW14R310J				300	J	
	APW14R010P				190	Pre-shape	
	APW14R310P				300	Pre-shape	
ASAHI SION blue ES	AHW14R017S	0.36 (0.014)	3	8.5	190	Straight	
	AHW14R317S				300	Straight	
	AHW14R017J				190	J	
	AHW14R317J				300	J	
	AHW14R017P				190	Pre-shape	
	AHW14R317P				300	Pre-shape	
Grand Slam	AG141002	0.36 (0.014)	4	4	180	Straight	
	AG141302				300	Straight	
	AG141002J				180	J	

Ordering Information & Lineup — CTO Guide Wires

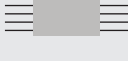Product	Catalog No.	Diameter (mm (inch))	Radiopaque Segment (cm)	Spring Coil Length (cm)	Length (cm)	Tip Shape	Clip Colour
Fielder XT-R	APW14R005S	0.26 / 0.36 (0.010 / 0.014)	16	16	190	Straight	
	APW14R305S				300	Straight	
Fielder XT-A	APW14R009S	0.26 / 0.36 (0.010 / 0.014)	16	16	190	Straight	
	APW14R309S				300	Straight	
ASAHI Gladius EX	AP14R025P	0.36 (0.014)	3	11	200	Pre-shape	
	AP14R325P				300		
ASAHI Gladius MG	AP14R023P	0.36 (0.014)	3	8.5	190	Pre-shape	
	AP14R323P				300		
ASAHI Miracle Neo 3	AH14R012P	0.36 (0.014)	10.5	10.5	190	Pre-shape	
	AH14R312P				300		
ASAHI Gaia Next 1	AH14R019P	0.36 / 0.27 0.014 / 0.011	15	15	190	Pre-shape	
	AH14R319P				300	Pre-shape	
ASAHI Gaia Next 2	AH14R020P	0.36 / 0.30 0.014 / 0.012	15	15	190	Pre-shape	
	AH14R320P				300	Pre-shape	
ASAHI Gaia Next 3	AH14R021P	0.36 / 0.30 0.014 / 0.012	15	15	190	Pre-shape	
	AH14R321P				300	Pre-shape	
ASAHI Conquest Pro / Confianza Pro	AGH143190	0.23 / 0.36 (0.009 / 0.014)	20	20	190	Straight	
	AGH143190P					Pre-shape	
	AGH143390				300	Straight	
	AGH143390P					Pre-shape	
ASAHI Conquest Pro 12 / Confianza Pro 12	AGH143191	0.23 / 0.36 (0.009 / 0.014)	20	20	190	Straight	
	AGH143191P					Pre-shape	
	AGH143391				300	Straight	
	AGH143391P					Pre-shape	
ASAHI Miracle 12 / MIRACLEbros 12	AG14M070	0.36 (0.014)	11	11	180	Straight	
ASAHI RG3	AHW10S302S	0.26 (0.010)	3	8	330	Straight	
ASAHI SUOH 03	AHW14R013P	0.36 (0.014)	3	19	190	Pre-shape	
	AHW14R313P				300	Pre-shape	
ASAHI EXTENSION 165cm	AG149001	0.36 (0.014)	—	—	165	—	

©2025 ASAHI INTECC CO., LTD. "ASAHI", "ASAHI SION", "Fielder", "Gladius" and "ULTIMATEbros" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.
Coated with SLIP-COAT® coating. Some materials/graphics were created in partnership with VP MED Group - www.vpmed.com

Manufactured by:

ASAHI INTECC
ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
Contact phone number: +81-561-86-9101
www.asahi-intecc.com

Contact:

ASAHI INTECC EUROPE B.V.
ASAHI INTECC Europe B.V.

Strawinskylaan 1461, 1077 XX Amsterdam, The Netherlands
Contact phone number: +31-20-794-0642
www.asahi-intecc.com