


Super Sheath™

Designed for Smooth Entry
and Reduced Trauma

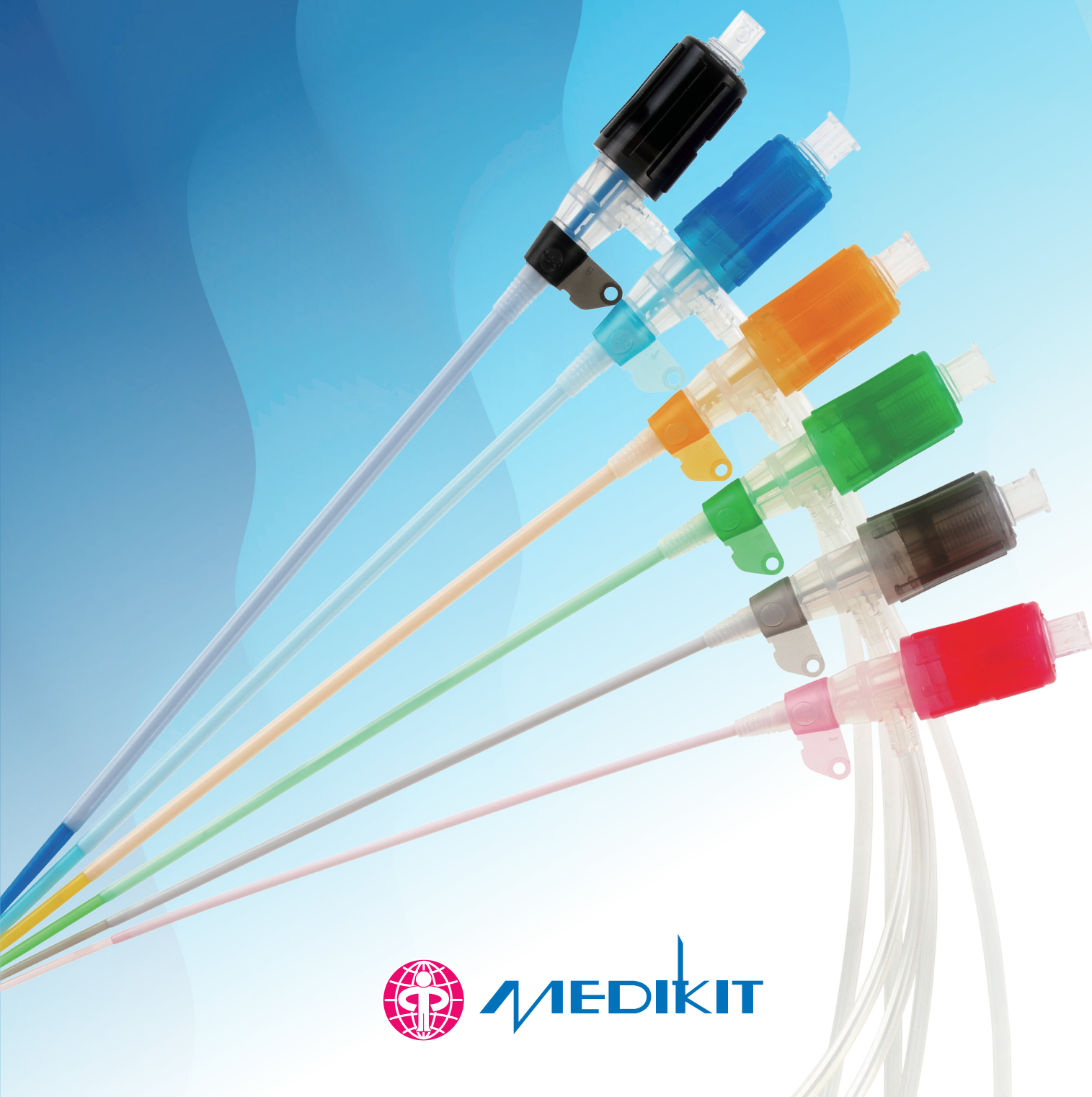


MEDIKIT

Super Sheath™

Designed for Smooth Entry and Reduced Trauma


Smooth Transitions

Funnel shaped dilator

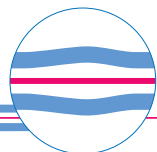

- ▶ The unique design of the dilator enables smooth transitions with a lower insertion force.
- ▶ Improved kink resistance and lubricity enhance pushability while maintaining device integrity.
- ▶ Long, gradually tapered tip reduces potential for vessel trauma.

Long, gradually tapered tip enables smooth device insertion

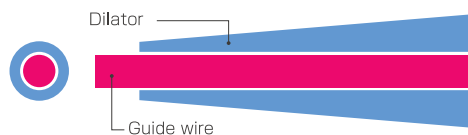


Excellent Catheter Trackability


- ▶ Silicone-coated valve and smooth inner surface promote device passage.
- ▶ Tricuspid valve design promotes hemostasis.

Ease of Use


Translucent hub

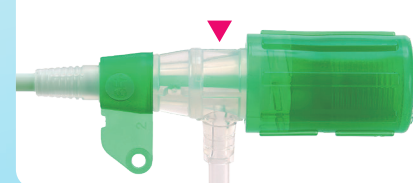

Rotating suture wing

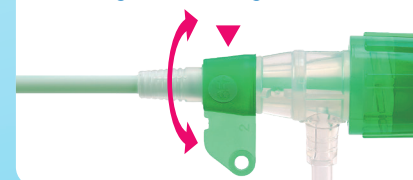


One-piece structure

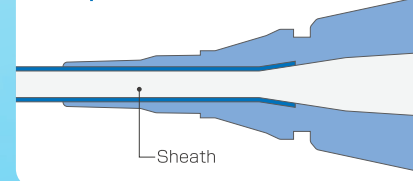


- ▶ Translucent hub allows visualization of the blood backflow.
- ▶ Rotating suture wing allows for sheath securement.
- ▶ Dilator/hub twist-lock for security during insertion.
- ▶ The one-piece structure of the sheath reduces device insertion/extraction resistance.

Various Sheaths for Various Procedures


Short (7cm,11cm) and Long (25cm) length Sheath

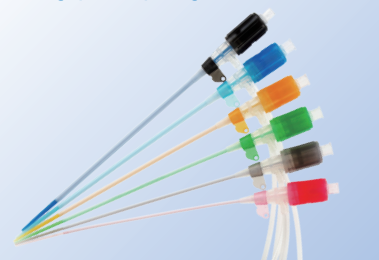

Larger French Sheath (10F, 11F, 12F, 14F)

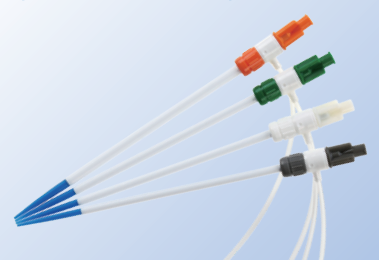


Sheath with a Radiopaque Marker

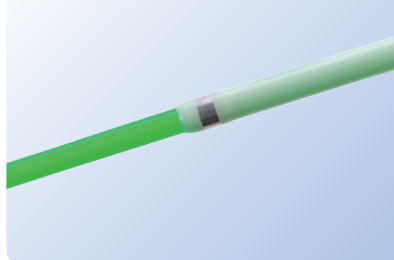

Super Sheath™

Components

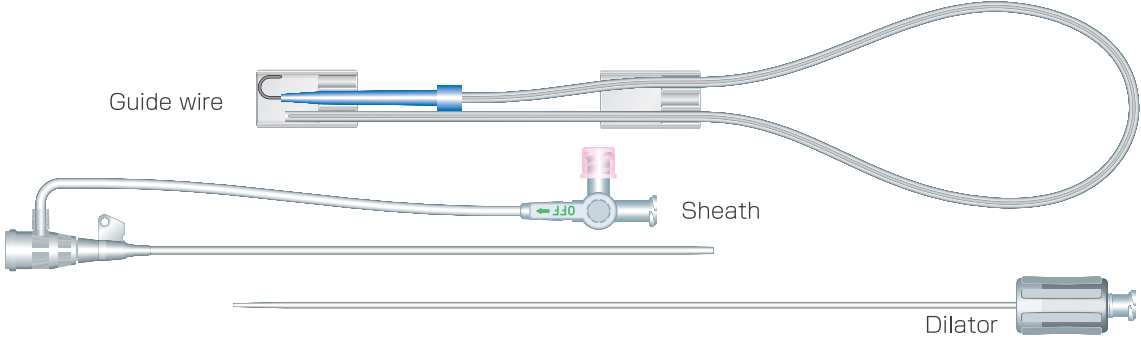

Standard Sheath

French Size	Sheath Length (cm)	Mini-Guide Wire(inch)
4F (red)	7, 11	0.035
5F (gray)	7, 11, 25	0.035, 0.038
6F (green)	7, 11, 25	0.035, 0.038
7F (orange)	7, 11, 25	0.035, 0.038
8F (blue)	11, 25	0.035, 0.038
9F (black)	11, 25	0.035, 0.038

Radiopaque Marker Sheath

Features	French Size	Sheath Length (cm)	Mini-Guide Wire (inch)
With Radiopaque Marker	4.0F (red)	11, 25	0.035
With Radiopaque Marker	5.0F (gray)	7, 11, 25	0.038
With Radiopaque Marker	6.0F (green)	7, 11, 25	0.038
With Radiopaque Marker	7.0F (orange)	7, 11, 25	0.038
With Radiopaque Marker	8.0F (blue)	7, 11, 25	0.038
With Radiopaque Marker	9.0F (black)	7, 11, 25	0.038

XL Sheath

Features	French Size	Sheath Length (cm)	Mini-Guide Wire (inch)
XL Sheath	10 (gray), 11 (white), 12 (green), 14 (orange)	11	0.035, 0.038
XL Sheath	10 (gray), 11 (white), 12 (green), 14 (orange)	25	0.038

- * All sheaths in the list above are packaged with dilators in a sterile pouch. One shelf box contains 10 sterile pouches.
- * 11cm sheaths may be purchased with or without Mini-Guide wires. 25cm sheath does not come with guide wires.
- * Sheaths with a radiopaque marker are supplied without a guide wire.



HEAD OFFICE: 1-13-2, Yushima, Bunkyo-ku, Tokyo 113-0034 JAPAN
<https://www.medikit.co.jp/english/contact/index.php>
MANUFACTURER: TOGO MEDIKIT CO., LTD. 17148-6, Hichiya, Hyuga-shi, Miyazaki, 883-0062, Japan

Please contact our local sales representative or agents to find out more details.
ASAHIINTECC USA, INC. Tel:1-949-756-8901, <http://www.asahi-inteccusa-medical.com/>
ASAHI INTECC EUROPE B.V. Tel:31-(0)20-794-0640 <https://www.asahi-intecc.eu/>